

TruKlear® – Rundum nur Vorteile

TruKlear® – Nothing but advantages



Die größten Vorteile auf einen Blick.

- | | |
|--|--|
| <p>1 Trapezförmiger, slotseitig breiter Verschlusschieber aus verfärbungs-freier Keramik für gute Rotationskontrolle.</p> | <p>4 Körper aus transluzenter Keramik mit abgerundeten Kanten für hohen Tragekomfort.</p> |
| <p>2 Abgerundete Slotteinläufe und vier Kontaktrippen im Slot reduzieren Klemm- und Kerbeffekte.</p> | <p>5 Metallfreier Haltemechanismus des Schiebers – für optimale Ästhetik, auch bei Nickelunverträglichkeit.</p> |
| <p>3 Vierwandslot ist besonders niedrig gestaltet, um Bögen mit einer Höhe von .025“ für eine gute Übertragung von Torque und Rotation eng zu führen.</p> | <p>6 Inverse Hakenbasis sorgt für optimale mechanische Retention und reduziert die Gefahr von Schmelzausrissen.</p> |

Most attractive advantages at a glance.

- | | |
|--|--|
| <p>1 Trapezoid, wide locking slider made of discoloration-proof ceramic for good rotation control.</p> | <p>4 Body made of translucent ceramic with rounded edges for excellent intraoral comfort.</p> |
| <p>2 Rounded slot inlets and four contact ribs in the slot reduce binding and notching.</p> | <p>5 Metal-free holding mechanism of the slider - for optimal aesthetics, even with nickel intolerance.</p> |
| <p>3 Four-wall slot is designed with a particularly low profile to precisely guide archwires with a height of .025“ ensuring effective transfer of torque and rotation.</p> | <p>6 Inverted hook base ensures perfect mechanical retention and reduces the risk of enamel fractures.</p> |


TruKlear® Brackets

TruKlear® brackets



Roth Prescription

Farbmarkierung | color-marks

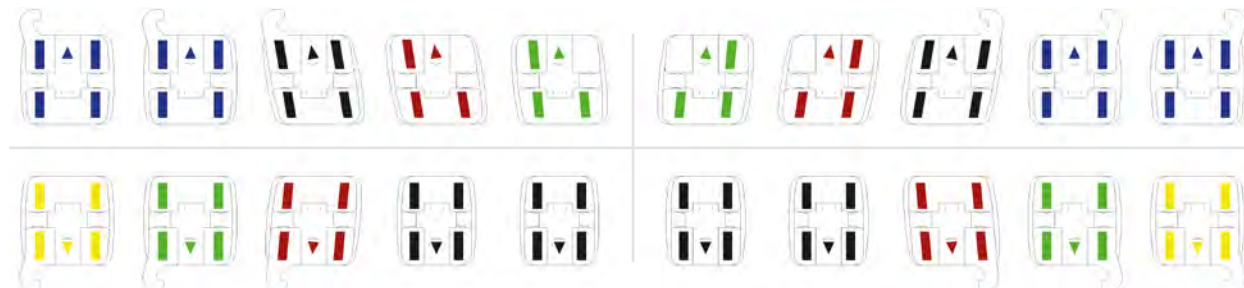

TruKlear® Brackets Roth Prescription

Maxillary					Slot .018"		Slot .022"	
Tooth	Torque	Angulation	In/Out	Rotation	Order No.		Order No.	
					Right	Left	Right	Left
1 Centrals	+12°	+5°	0,9	–	C713-0103	C712-0103	C713-0101	C712-0101
2 Laterals	+8°	+9°	1,1	–	C713-0203	C712-0203	C713-0201	C712-0201
3 Cuspids + hook	-2°	+11°	0,7	4° mesial	C713-0313	C712-0313	C713-0311	C712-0311
4 Bicuspid	-7°	0°	0,7	2° distal	C713-0503	C712-0503	C713-0501	C712-0501
4 Bicuspid + hook	-7°	0°	0,7	2° distal	C713-0513	C712-0513	C713-0511	C712-0511
5 Bicuspid	-7°	0°	0,7	2° distal	C713-0503	C712-0503	C713-0501	C712-0501
5 Bicuspid + hook	-7°	0°	0,7	2° distal	C713-0513	C712-0513	C713-0511	C712-0511

Mandibular					Slot .018"		Slot .022"	
Tooth	Torque	Angulation	In/Out	Rotation	Order No.		Order No.	
					Right	Left	Right	Left
1 Centrals	0°	0°	1,4	–	C712-1203	C712-1203	C712-1201	C712-1201
2 Laterals	0°	0°	1,4	–	C712-1203	C712-1203	C712-1201	C712-1201
3 Cuspids + hook	-11°	+7°	1,2	2° mesial	C713-1413	C712-1413	C713-1411	C712-1411
4 Bicuspid + hook	-17°	0°	1,2	4° distal	C713-1513	C712-1513	C713-1511	C712-1511
5 Bicuspid + hook	-22°	0°	1,2	4° distal	C713-1613	C712-1613	C713-1611	C712-1611

Fallsätze | Cases

Variation	Slot .018"			Slot .022"		
Cases	1	3	6	1	3	6
 3 Max.	C706-2000	C706-2001	C706-2002	C706-2006	C706-2007	C706-2008
 3-5 Max.	C706-2003	C706-2004	C706-2005	C706-2009	C706-2010	C706-2011
 3-5 Max. + Mand.	C706-2027	C706-2028	–	C706-2033	C706-2034	–